

FENTANYL TEST STRIP USER GUIDE

1. PREPARE

INJECTING

- Prepare your shot
- After preparing, add 10 drops* of water to your cooker and stir well

PILLS/POWDER

- Chop/crush drugs into fine powder & mix well
- Add water* to a baggie, cup, or water bottle cap
- Mix **50 mg** of your drugs into the water
- Stir well

*About a teaspoon

2. OPEN TEST STRIP

- Hold the **solid blue** end of the test strip
- Dip into water up to the blue line for **15 seconds** then take it out
- Wait for **3 minutes** to check your results



1

10 mg = enough to cover the hair of Abraham Lincoln's head on a penny


10mg for MDMA & Meth

50mg for all other drugs


3

3. READ RESULTS

No water past this line

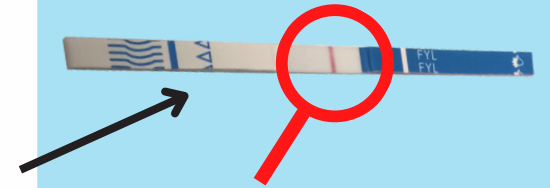


Positive for fentanyl:
ONE RED line

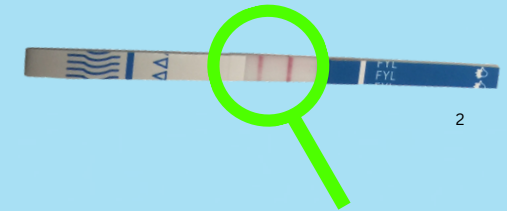

Negative for fentanyl:
TWO RED lines

The second line can be faint.
This still means negative.

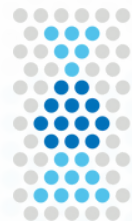

MATTERS[®]
Medication for Addiction Treatment & Electronic Referrals

1. NYC Health. (n.d.). How to test your drugs using fentanyl test strips <https://www1.nyc.gov/assets/doh/downloads/pdf/basas/fentanyl-test-strips-brochure.pdf>

2. raveready. (2019). Fentanyl Test Strips Kit. <https://www.raveready.com/dietary-health/3828-fentanyl-test-strips-kit.html>; 3. DanceSafe. (2022). How to Test Your Drugs for Fentanyl. [Infographic] dancesafe.org